




## Product Overview

Safe n Simple focuses on product quality, patient safety, and cost containment. With the introduction of our 100% collagen wound care dressings we are proud to launch the Simpurity wound care product line. Simpurity wound care products follow evidence based medicine, facilitate moist wound healing and improve patient outcomes while lowering costs.

Simpurity Collagen wound dressings consist of lyophilized medical grade bovine collagen presenting as a clean, white, single layered sponge. Simpurity Collagen undergoes an exclusive purification process and is highly biocompatible as a result. And, like human connective tissue, Simpurity Collagen wound dressings are rich in Type One Collagen with a unique porous structure (photo at right) that provides excellent absorption and scaffolding for cell proliferation and migration.

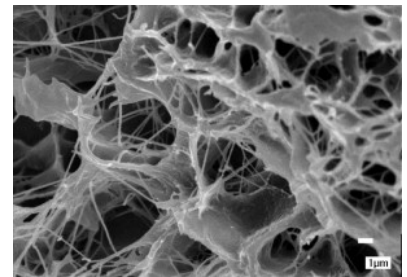

Simpurity Collagen provides substrate collagen to divert matrix metalloproteinases (MMP) digestion from healing tissues and reduces inflammation that may slow or prevent healing. Simpurity Collagen fragments provide ingredients for tissue regeneration and the 3D scaffold structure allows fibroblasts to migrate across the wound bed facilitating granulation, angiogenesis and tissue repair. In acute wounds, Simpurity Collagen rapidly activates platelets and coagulating factors to induce hemostasis.



## Product Features

- Low antigenic collagen will minimize an immune reaction
- Glutaraldehyde-free crosslinking process reduces immunogenicity
- Atelocollagen enhances low inflammatory potential
- Soft and highly flexible for easy application and superior conformity
- Creates a soft conformable gel on contact with wound fluid
- Highly biocompatible porous scaffold encourages healing
- Crosslinking process increases mechanical strength and moderates enzymatic digestion
- High exudate absorption and retention provides moist wound healing environment
- May be used as a hemostatic dressing
- 100%, non-bleached, native undigested Bovine Collagen

## Indications

- Partial and full-thickness wounds
- Tunneled / undermined wounds or surgical wounds
- Low-to-moderately exuding chronic wounds including
  - Diabetic foot ulcers
  - Venous leg ulcers
  - Pressure sores
- Management of light bleeding in acute or surgical wounds
- Post debridement for hemostasis

## Product Details

Product	Item Number	Units/Box	Boxes/Case	HCPC
2" x 2"	SNS50002	10	14	A6021



Clarkston, MI 48348 | 844-767-6334 | [www.SnS-Medical.com](http://www.SnS-Medical.com)

WCC1-030116